



# FTIR engine (FT-NIR spectrometer) Measurement Examples

## CONTENTS

Transmission Measurements	3
- Principle of Transmission Method	P.4
- Examples of Transmission Measurement (Water, Ethanol)	P.5
- Concentration Analysis of Alcoholic Beverages	P.6
- Concentration Analysis of Liquid Chemicals	P.7
Diffuse Reflection Measurements	9
- Principle of Diffuse Reflection Method	P.10
- Moisture Content Measurements	P.11
- Identification of Powder Samples	P.12
- Identification of Additives in Various Plastics	P.13
- Analysis of Concrete Degradation	P.14
- Film Thickness Measurements	P.15
- Pharmaceutical Analysis	P.16

# Near infrared spectroscopic analysis using an FTIR engine


Molecules each have unique vibrations, thereby absorbing near infrared light of a specific wavelength. Near infrared spectroscopic analysis utilizes this characteristic to analyze components contained in substances. In the near infrared region of 1100 nm to 2500 nm, many substances have unique absorption spectra, and these are applied to qualitative and quantitative analysis in various fields.

The FTIR engine is a compact Fourier transform infrared spectrometer developed for near infrared spectroscopic analysis. A Michelson optical interferometer and control circuit are built into a palm-sized enclosure. Spectrum and absorbance can be measured by connecting a PC via USB. It can be applied to real-time measurement performed on site without bringing the measurement sample into the analysis room as well as continuous monitoring.

Two types of measurement methods of "transmission measurement" and "reflection measurement" are used for near infrared spectroscopic analysis using the FTIR engine.


## Optical system

The Michelson interferometer is used for the FTIR engine. The FTIR's optical interferometer is composed of a light input section, beam splitter, photodetector, and MEMS chip. The MEMS chip has a movable mirror ( $\phi 3$  mm) that uses MEMS (micro-electro-mechanical systems) technology and a fixed mirror. The photodetector (InGaAs PIN photodiode) acquires light intensity signals that varies depending on the movable mirror position. The optical spectrum is obtained by taking the Fourier transform of this light intensity signal. The built-in semiconductor laser (VCSEL: vertical cavity surface emitting laser) and photodetector (Si photodiode) for monitoring the movable mirror position allow spectrum measurement with high wavelength accuracy.

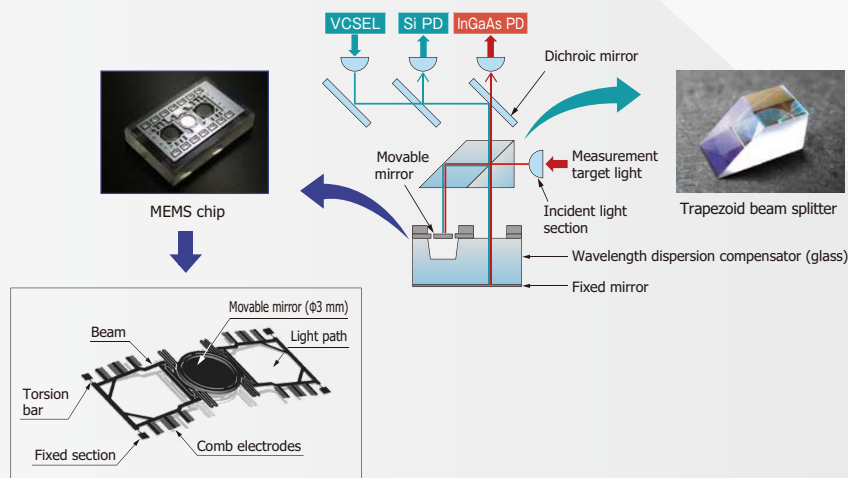


### MEMS actuator

The movable mirror uses an electrostatically driven MEMS actuator. The mirror is moved up and down in parallel by applying voltage to comb electrodes.

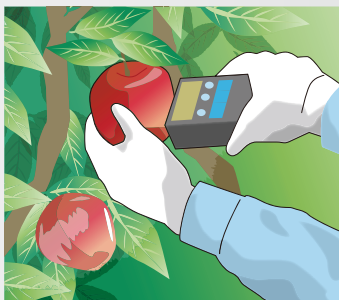


Pre-harvest inspection of agricultural products

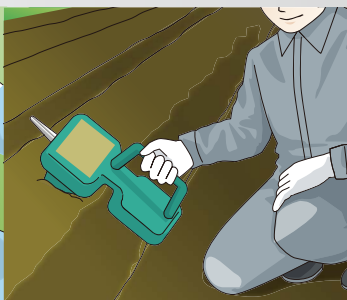

Soil analysis



Material acceptance inspection

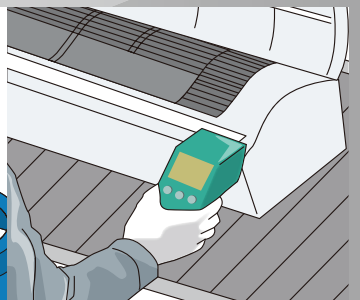

Plastic sorting



# Transmission Measurements

# Principle of Transmission Method

Near-infrared spectroscopy uses an analytical method that relates band intensity of a spectrum (absorbance) to sample concentration.

Absorbance is defined by the equation  $A_\lambda = -\log_{10} (I_1/I_0)$ . And common logarithms of a ratio between incident light level  $I_0$  (reference measurement) and transmitted light level  $I_1$  (sample measurement) can be taken for analysis.

While transmittance decays exponentially with increasing optical path length, absorbance expressed in logarithmic terms changes proportionally with optical path length.

For example, if transmittance is 0.1 (absorbance = 1) and the thickness of the object is doubled, transmittance will be 0.01, while absorbance will be doubled (absorbance = 2).

## Equation of absorbance

$$A = -\log_{10} \left( \frac{I_1}{I_0} \right) = cI\varepsilon$$

A: absorbance


c: sample concentration

I: optical path length of quartz cell

$\varepsilon$ : molar absorption coefficient


$I_1$ : transmitted light level

$I_0$ : incident light level

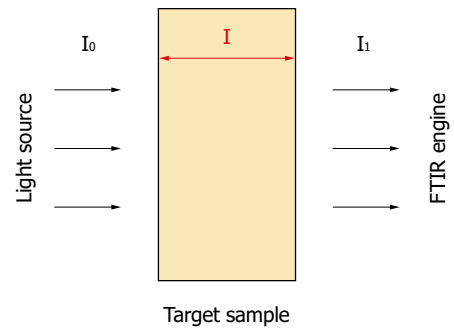

KACCC1126EA

## Absorption spectrum measurement system

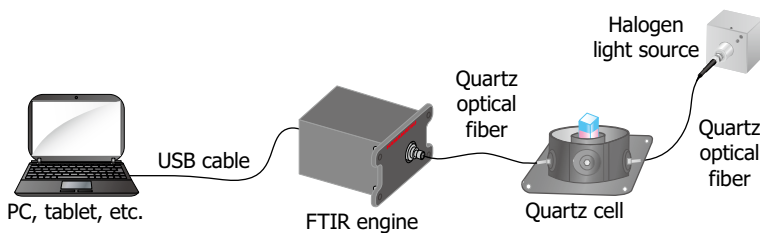


KACCC1127EA

# Examples of Transmission Measurements (Water, Ethanol)

## Water absorption spectrum


The water absorption spectrum has a peak of each 1450 nm (OH groups) and 1930 nm (H<sub>2</sub>O groups) band. The absorbance changes according to the optical path length of the quartz cell. Absorbance is measured well up to about 2.0.

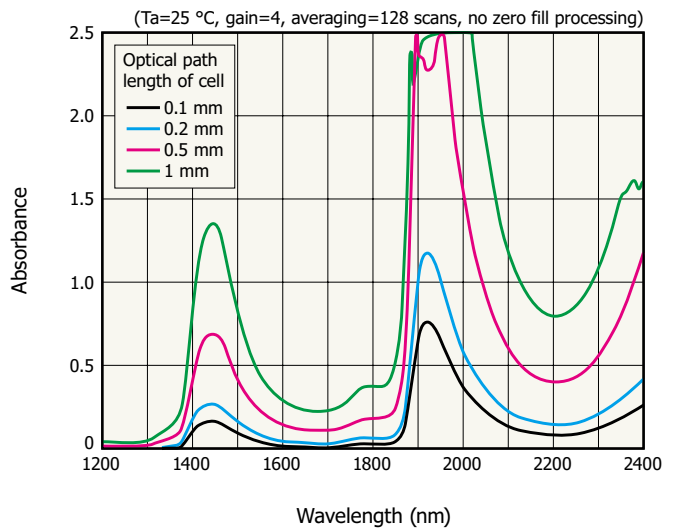

KACCB0614EA

## Ethanol absorption spectrum

The ethanol absorption spectrum has a peak after 2200 nm. Zero-fill processing is performed before the Fourier transform, and data points are interpolated to display a smooth spectrum.

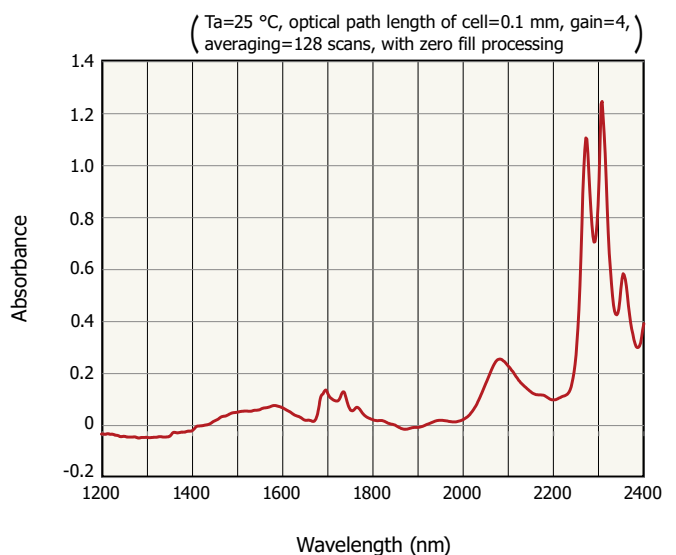


KACCB0693EA

# Concentration Analysis of Alcoholic Beverages

## Comparison of absorbance of alcoholic beverages and estimation of alcohol concentration


Following figure shows the near-infrared absorption spectra of beer, sake, brandy, ethanol, and water. There is absorption by the OH group of water (1450 nm band, 1900 nm band) and by the CH group of alcoholic beverages (2100 to 2500 nm). With transmission measurement results, we were able to obtain characteristic spectra in the absorption bands of water and alcoholic beverages.

In addition, with the results of estimating the alcoholic concentration from absorbance in the 2300 nm band, we confirmed that the estimated values and numerical values of components contained in the beverage matched, and that high accuracy measurement is possible.

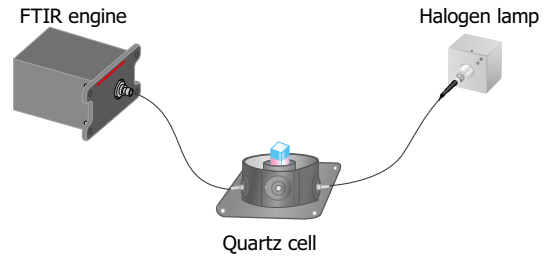

KACCC1128EA


### ■ Absorption spectra of alcohols (Measurement example)

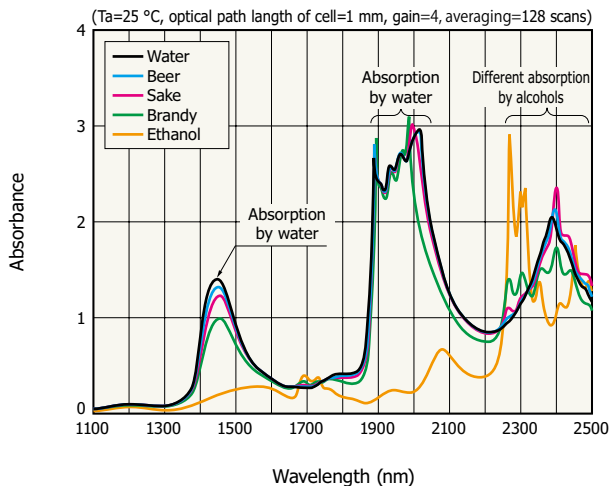

KACCB0622EB

Calibration curve based on ethanol measurement results

### ■ Alcohol concentrations (measurement example of 2300 nm band)

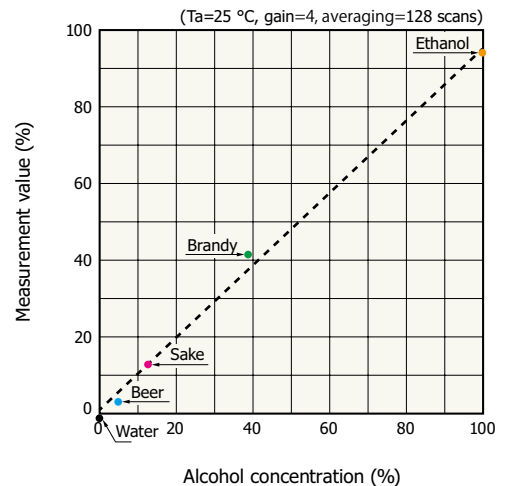


KACCB0623EB

Measurement cooperated by: Hamamatsu Central Research Laboratory

# Concentration Analysis of Liquid Chemicals

## Accurate quantitative measurements of various liquid chemical concentrations

Quantitative analysis can be done by measuring samples of liquids with different concentrations. Figures on page 7 and page 8 show the measurement examples of nitric acid aqueous solutions ( $\text{NO}_3^-$ ) with different concentration.

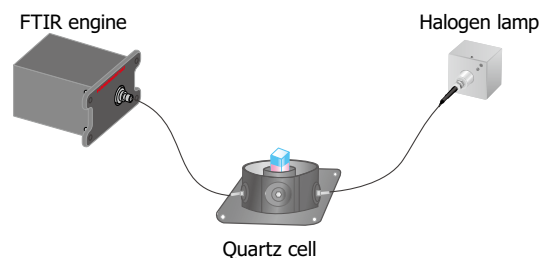


KACCC1128EA

### ■ Measurement samples

Measured using 6 samples each of 3 types (18 total) of nitric acid aqueous solutions with different concentrations

- ① High concentration ( $\text{NO}_3^-$ : 0%, 2%, 4%, 6%, 8%, 10%)
- ② Medium concentration ( $\text{NO}_3^-$ : 0%, 0.2%, 0.4%, 0.6%, 0.8%, 1%)
- ③ Low concentration ( $\text{NO}_3^-$ : 0%, 0.02%, 0.04%, 0.06%, 0.08%, 0.1%)


### ■ Absorption spectra of nitric acid in aqueous solution (High concentration)

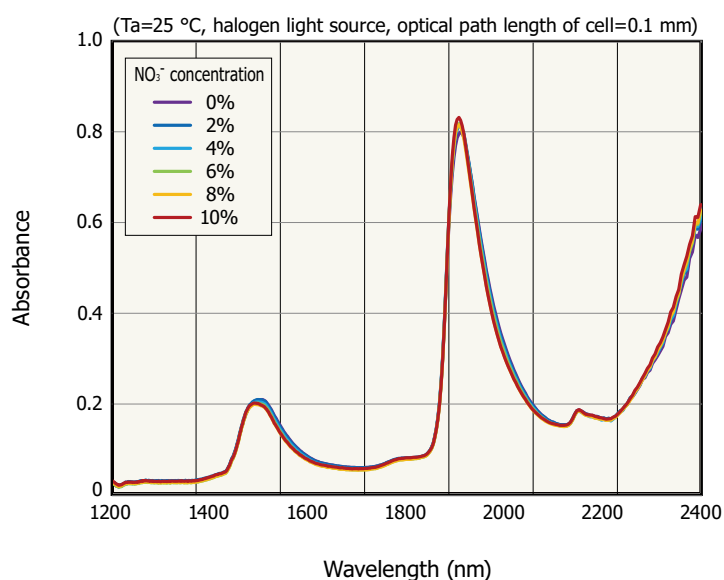

KACCB0700EA

Measurement cooperated by: Hamamatsu Central Research Laboratory

Following figures show regression coefficient and calibration curve of different concentration of nitric acid aqueous solution ( $\text{NO}_3^-$ ). In the high and medium concentration measurement samples, there is a large absorption change in the 2000 nm band (OH group of water), and there is a small change in the 2200 nm band (OH group and NH group of nitrogen). While for low concentration measurement samples, there are little change of absorption peaks at particular wavelength. A calibration curve is formed by using the data from the entire wavelength range. When the concentration is lower, quantification accuracy will also be lower, but the level is high enough to do concentration measurement.

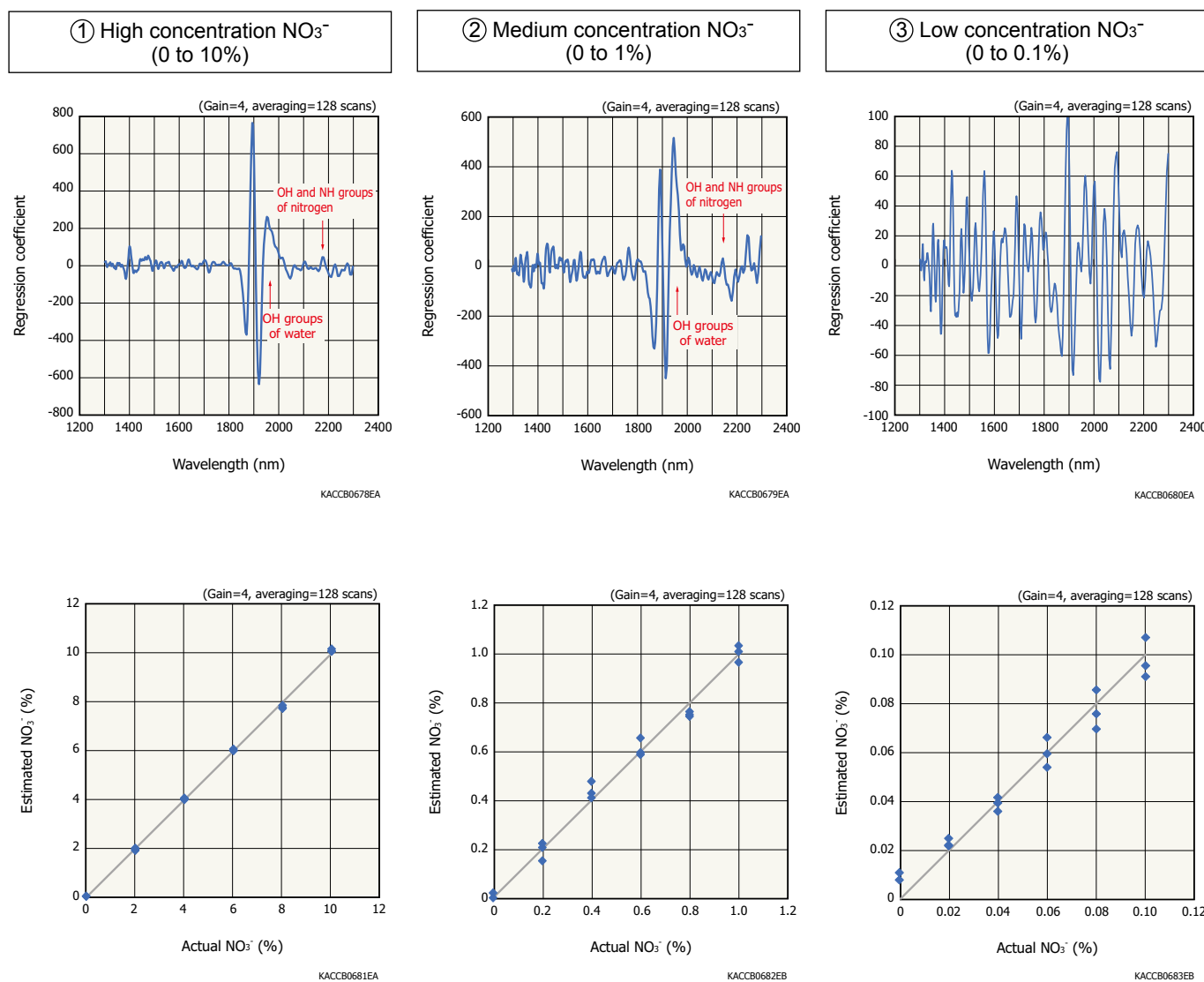

Measurement cooperated by: Hamamatsu Central Research Laboratory





# Diffuse Reflection Measurements

# Principle of Diffuse Reflection Method


While a portion of light irradiated onto a sample is regular reflected by surface of particles, the rest penetrates into the sample. The light is repeatedly diffused through refractive transmission, light scattering, surface reflection inside the sample, and some light emitted out from the sample's surface is measured. The diffuse reflection spectrum is similar to the absorption spectrum because the light is repeatedly transmitted through the interior of the sample during the light diffusion process. In diffuse reflection measurement, common logarithms of a ratio between incident light level  $I_0$  (reference measurement) and transmitted light level  $I_1$  (sample measurement) can be taken for analysis.

## Equation of absorbance

$$\frac{K}{S} = \frac{(1 - R)^2}{2R} = \cosh [\log_{10} (1/R)] - 1 \cong \log_{10} \left( \frac{I_1}{I_0} \right)$$

K/S: Kubelka-Munk  
 S: scattering coefficient  
 K: absorption coefficient

R: Reflectance =  $I_1 / I_0$   
 $I_1$ : transmitted light level  
 $I_0$ : incident light level

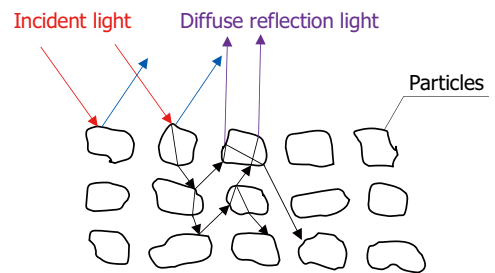


KACCC1129EB

## Related products

### Diffuse reflection light source L16462-01



Diffuse reflection light source integrates lamps and bundled optical fibers. Light from lamps is irradiated to a sample, and is diffused and reflected. The light which is guided into the optical fibers can be used for spectroscopic analysis by connecting with a near-infrared spectrometer. By arranging multiple lamps and optical fibers close to each other, it can effectively detect the weak diffused light emitted from the sample.

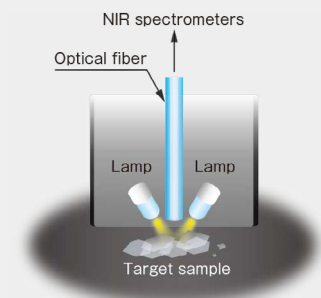

KACCC1131EA


[Structure of diffuse reflection light source]

# Moisture Content Measurements

## Measures absorption at 1450 nm and 1930 nm


We prepared a dry cloth and a wet cloth soaked in water, then compared their absorption spectra with the FTIR engine.

The wet cloth has a stronger absorption spectrum than the dry cloth in the 1450 nm band (OH group) and 1930 nm band (H<sub>2</sub>O group), which are the water absorption bands.

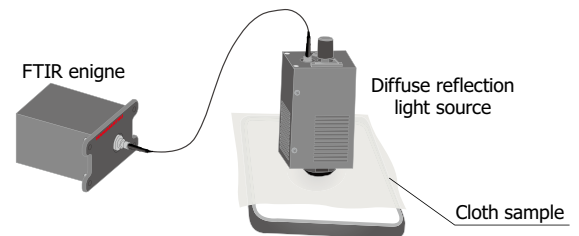

KACCC1132EA

### ■ Comparison of absorbance spectra of wet cloth and dry cloth samples

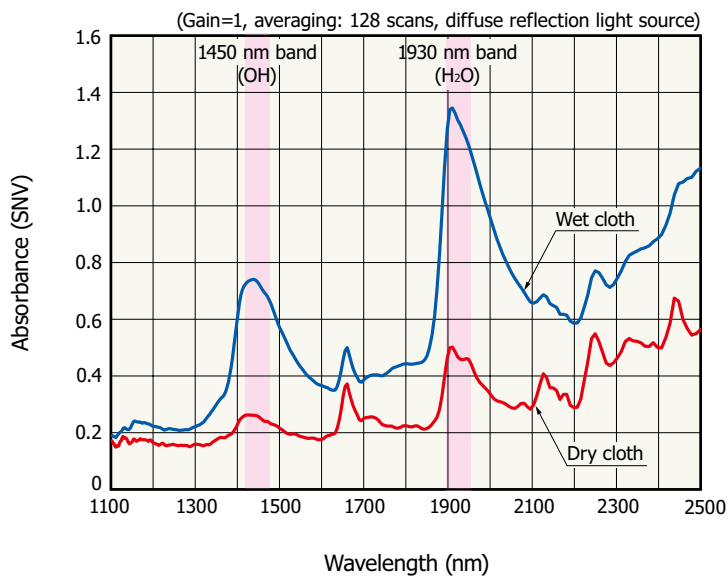

KACCB0697EA


SNV: Standard Normal Variate

# Identification of Powder Samples

## Comparison of absorbance of sugar


Sugars such as monosaccharides (glucose and fructose) and disaccharides (sucrose) can be found in food, defining its nutritional and flavor content.

The near infrared spectroscopy is applied for composition measurement of saccharide content in fruits. The FTIR engine is capable of accurately measuring the most minute peak patterns, resulting in spectra that are similar to those produced using a benchtop spectrometer.

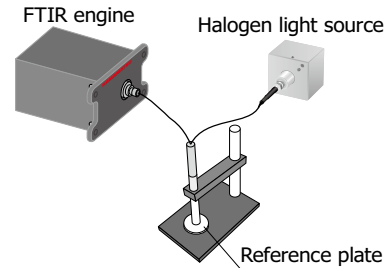

KACCC1133EA

## ■ Comparison of absorbance spectra of powder sugar samples

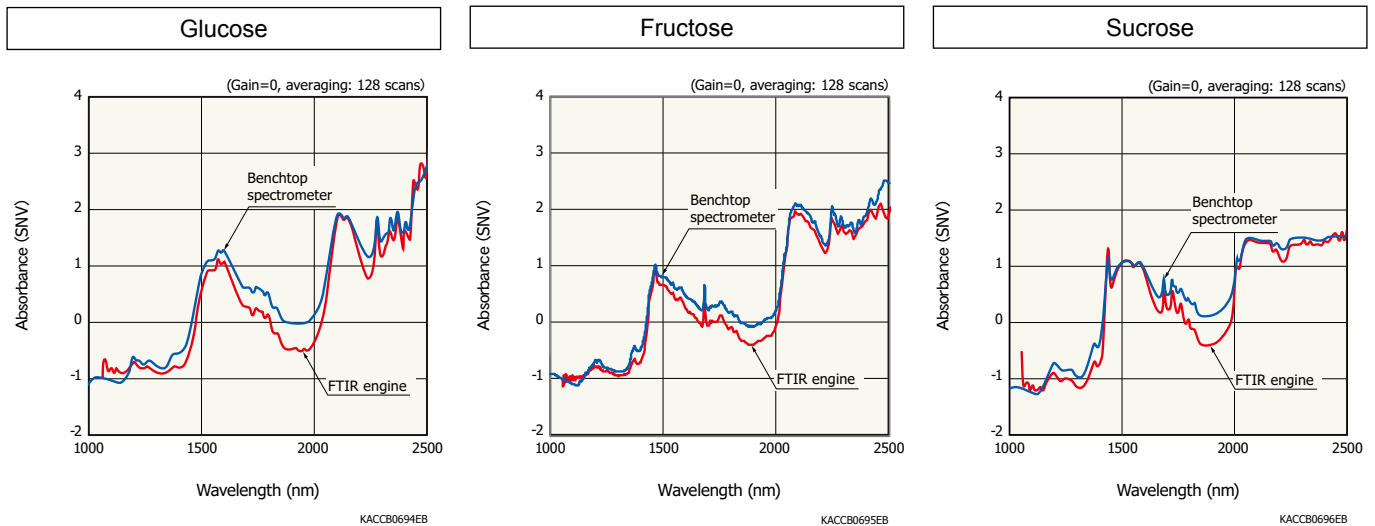


Measurement cooperated by: Hamamatsu Central Research Laboratory

# Identification of Additives in Various Plastics

## Classification of brominated flame retardants


Various types of plastics are used in consumer electronics. For the measurement of white plastic (PP: polypropylene), differences in spectra in a region from 2000 to 2500 nm make separation of brominated flame retardants used in their plastics possible.

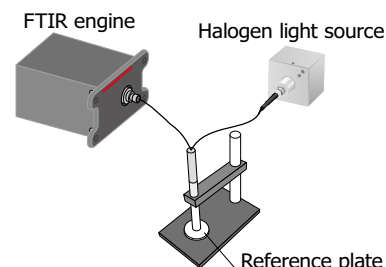

KACCC1133EA

### ■ Contents measurement of white polypropylene


The right figure shows absorbance spectra of additives such as DBDE (decabromodiphenyl ether), TBBA (tetrabromobisphenol A) and Talc. From the absorption peak at around 1400 nm, the graph shown in red is possible to infer a talc-specific component. Although TBBA and DBDE have relatively similar spectra, they can be distinguished from each other by the difference in the absorption peak at around 2100 nm.

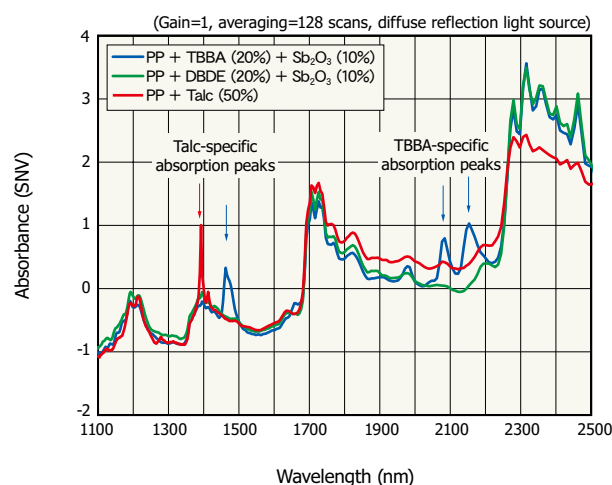

KACCB0698EA

### ■ Estimation of concentration level

The right figure shows absorbance spectra of white polypropylene integrating different concentration additives of TBBA and  $Sb_2O_3$  (antimony trioxide). The concentration of TBBA additives can be estimated from the difference in absorbance intensities at 1700 nm, 2200 nm and 2300 nm.

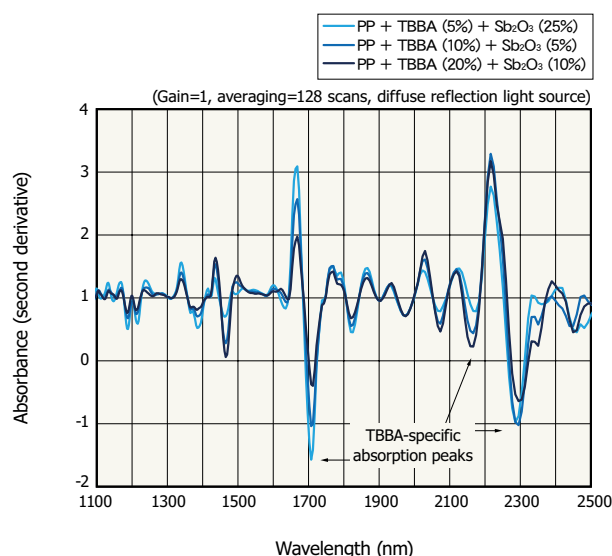

KACCB0699EA

Measurement cooperated by: Hirofumi Kawazumi  
Research Fellow,  
[Institute of Systems, Information Technologies and Nanotechnologies (ISIT)]

# Analysis of Concrete Degradation

## Estimating natural deterioration of concrete blocks


Following shows measurement results of concrete blocks installed outdoors for one and fifty years. Estimating natural deterioration of concrete blocks is possible by measuring progression of hydration and calcium (Ca) dissolution.

It also shows increase of the hydration (OH, H<sub>2</sub>O) around 1430 nm and 1930 nm, and CSH around 2210 nm.



Material age	50 years	18 years	10 years	1 year
Installation environment	Outdoors	Under a roof	Outdoors	Outdoors

### ■ Measurement of concrete blocks left outdoors (1 year vs. 50 years)

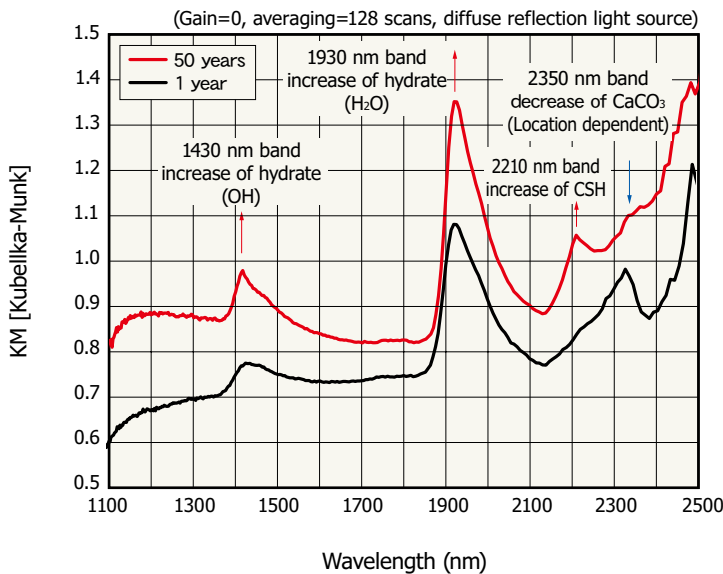

Measurement cooperated by: Emeritus Professor Satoru Nakashima (Osaka Univ.)



# Film Thickness Measurements

## Measuring film thickness by spectral interferometry


When white light enters a thin film sample, multiple reflections occur inside the film. These multiple-reflection light waves boost or weaken each other along with their phase difference. The phase difference of each multiple-reflection light is determined by the light wavelength and optical path length (distance that light moves back and forth in the thin film multiplied by the film refractive index).

This phase difference allows the spectrum reflected from or transmitted through the sample to produce a unique spectrum that depends on the film thickness. Spectral interferometry is a technique for measuring film thickness by analyzing that particular spectrum.

Following figure shows the result of film thickness measurement example of 2-layer plastic film and glass plate measured by the FTIR engine.

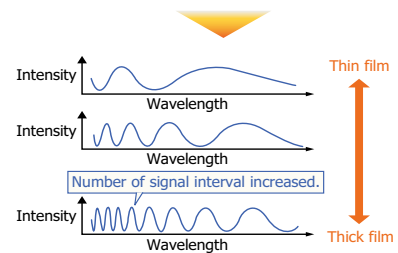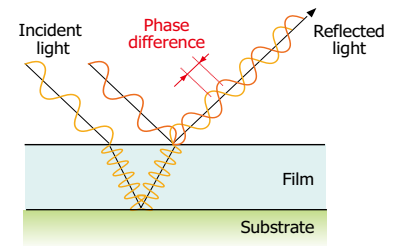

### ■ Film thickness measurement (2-layer plastic film, glass plate)

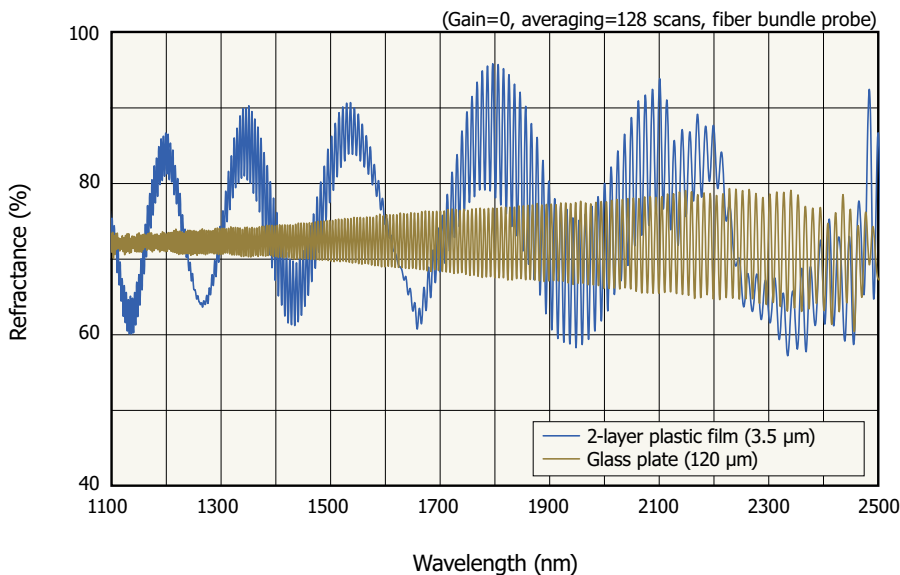

KACCB0702EA


# Pharmaceutical Analysis

## Absorbance measurement of medicines

Following shows spectra of five types of medicines including stimulants (Caffeine, Theophylline) and analgesics (Loxoprofen, Paracetamol, Aspirin).


Measurement results of Theophylline component using an FTIR engine show peaks in the specific spectrum is assigned and consistent with a published paper <sup>1)</sup>.

1) K Ikegaya, et al., Nippon Shokuhin Kogyo Gakkaishi, 34(4), p255(1987)

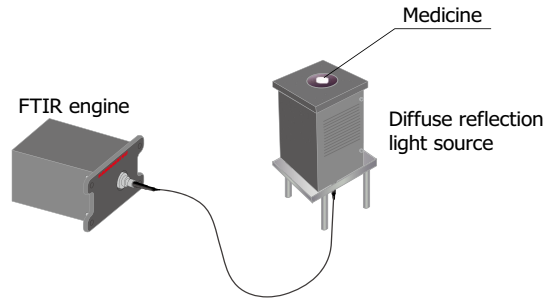

KACCC1134EA

## ■ Measurement of medicines

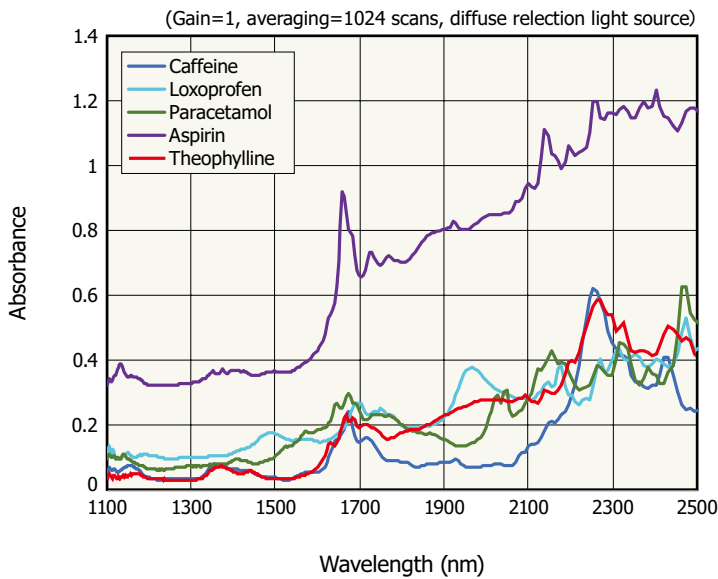


KACCB704EA

Measurement cooperated by: Hamamatsu Central Research Laboratory



## Quantitative analysis of Theophylline Anhydrate

Tablets containing theophylline, caffeine, and additives were analyzed using Hamamatsu's diffuse reflection light source L16462-01 and FTIR engine C15511-01. Using the same measurement samples (6 tablets with different theophylline concentrations), only theophylline was quantitatively measured. Theophylline content is predicted by making a calibration curve from previously acquired training data\*<sup>1</sup> and doing PLS\*<sup>2</sup> regression analysis. The coefficient of determination on the calibration curve is  $R^2=0.997$ , showing high accuracy measurements are possible. Even in conditions in which the tablets contain caffeine and theophylline, which have similar spectra in the near-infrared region, it is possible to quantitatively analyze only theophylline.

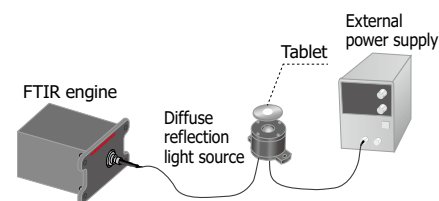

KACCC1153EC

\*1: When quantifying a substance, this is data that expresses the change with respect to concentration by measuring absorbance of the standard substance in advance.

\*2: Partially Square Regression


### Measurement samples

Sample No.	t0	t10	t20	t30	t40	t50
Theophylline	0%	10%	20%	30%	40%	50%
Caffeine	20%	20%	20%	20%	20%	20%
Additives (Starch, Lactose)	80%	70%	60%	50%	40%	30%

### Measurement conditions

- Program: Uses Python's scikit-learn module
- Samples with 0% content (t0) were excluded from the analysis to determine the MAPE (mean absolute percentage error)
- Number of PLS principal ingredients : 3


### Absorption spectra of the components included in the tablet

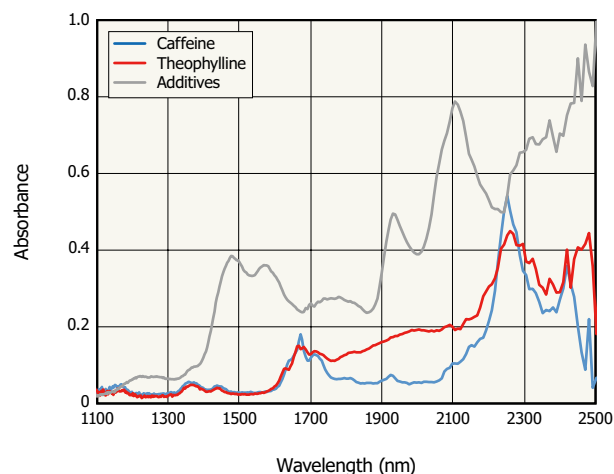

KACCB0742EA


\* Caffeine and theophylline are close in spectrum in the near-infrared region 1100 nm to 1700 nm.

### Absorbance spectra of the tablets

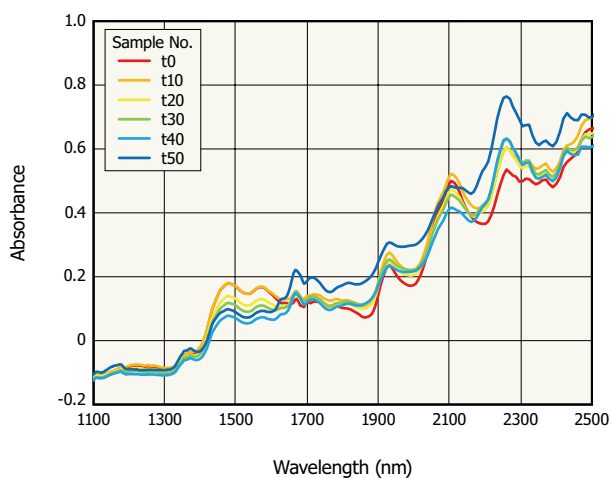

KACCB0734EA

### Calibration curve of Theophylline

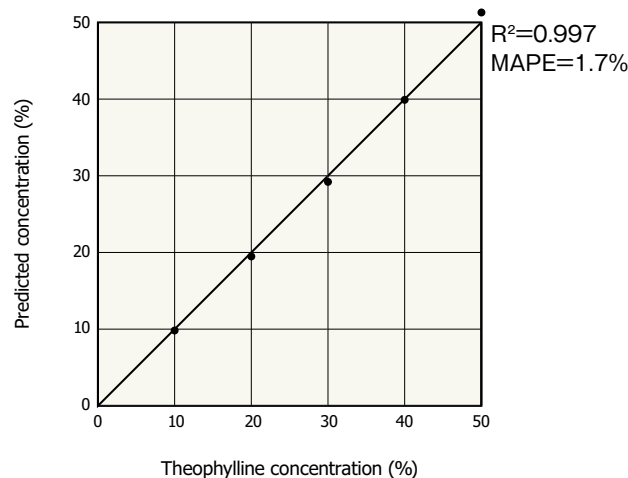


KACCB0664EA

Measurement cooperated by: Hamamatsu Central Research Laboratory

# "Zero Fill Processing" Enables Spectrum Detailed Analysis


Zero fill processing is a technique of processing the obtained spectrum data to finer detail. It is a process of adding zeros to each end of an optical interference signal before a Fourier transform is applied. This allows for interpolation between points that are plotted after the Fourier transform. The red graph using the zero fill processing achieves similar data compared with the reference data.

## ■ Ethanol absorption spectrum

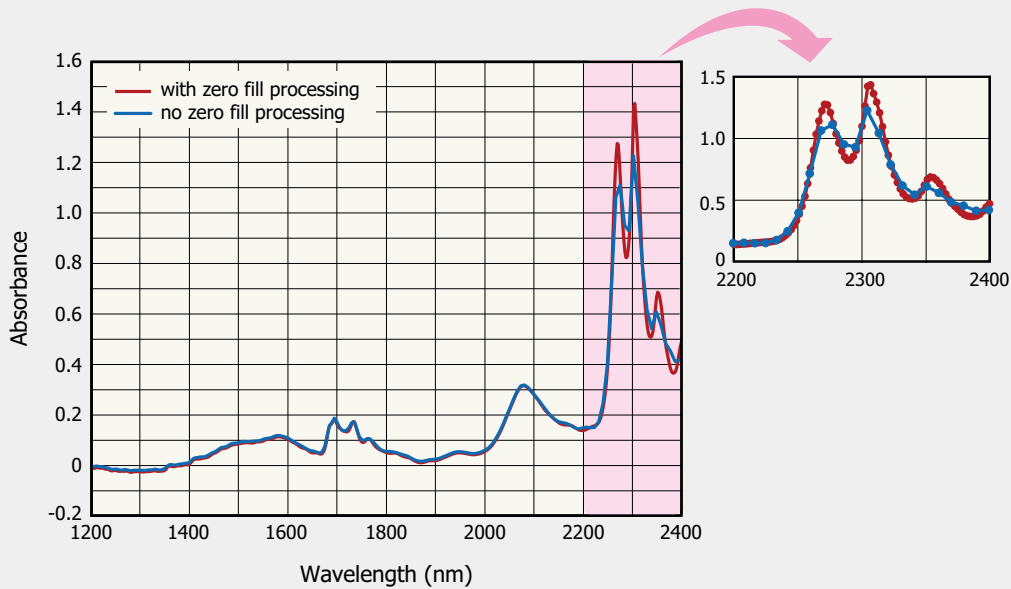

KACCB0706EA


## □ Reference: literature data

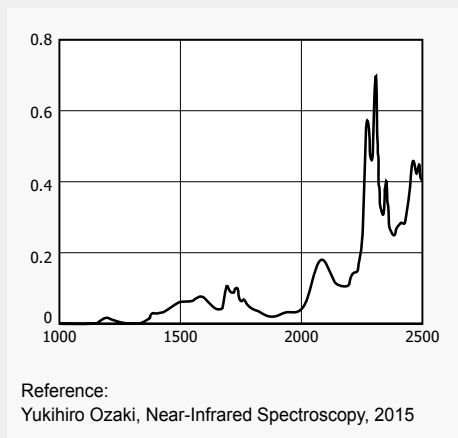

# "Reference Measurement" Realizing Stability under Continuous Operation


In the FTIR engine, the internal temperature heats up 1 hour after the start of operation and then stabilizes. Output variation is suppressed by performing sample measurement and reference measurement alternately.

## ■ FTIR internal temperature and output variation during reference measurement

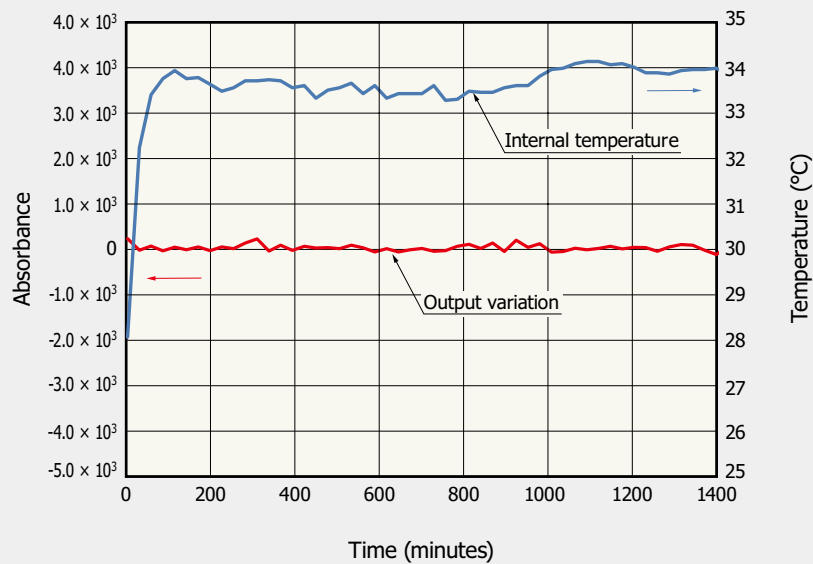

KACCB0674EA

## ■ Connection example of measurement system

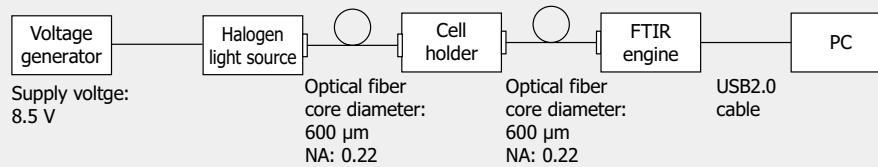

KACCC1122EA

## ■ Technical notes

FTIR engine (FT-NIR spectrometer)

[https://www.hamamatsu.com/resources/pdf/ssd/ftir\\_engine\\_kacc9012e.pdf](https://www.hamamatsu.com/resources/pdf/ssd/ftir_engine_kacc9012e.pdf)

Information described in this material is current as of September 2024.

Product specifications are subject to change without prior notice due to improvements or other reasons. This document has been carefully prepared and the information contained is believed to be accurate. In rare cases, however, there may be inaccuracies such as text errors. Before using these products, always contact us for the delivery specification sheet to check the latest specifications.

The product warranty is valid for one year after delivery and is limited to product repair or replacement for defects discovered and reported to us within that one year period. However, even if within the warranty period we accept absolutely no liability for any loss caused by natural disasters or improper product use. Copying or reprinting the contents described in this material in whole or in part is prohibited without our prior permission.

# HAMAMATSU

[www.hamamatsu.com](http://www.hamamatsu.com)

HAMAMATSU PHOTONICS K.K., Solid State Division

1126-1 Ichino-cho, Chuo-ku, Hamamatsu City, 435-8558 Japan, Telephone: (81)53-434-3311, Fax: (81)53-434-5184

U.S.A.: HAMAMATSU CORPORATION: 360 Foothill Road, Bridgewater, NJ 08807, U.S.A., Telephone: (1)908-231-0960, Fax: (1)908-231-1218

Germany: HAMAMATSU PHOTONICS DEUTSCHLAND GMBH: Arzbergerstr. 10, 82211 Herrsching am Ammersee, Germany, Telephone: (49)8152-375-0, Fax: (49)8152-265-8 E-mail: [info@hamamatsu.de](mailto:info@hamamatsu.de)

France: HAMAMATSU PHOTONICS FRANCE S.A.R.L.: 19 Rue du Saule Trapu, Parc du Moulin de Massy, 91882 Massy Cedex, France, Telephone: (33)1 69 53 71 00, Fax: (33)1 69 53 71 10 E-mail: [infos@hamamatsu.fr](mailto:infos@hamamatsu.fr)

United Kingdom: HAMAMATSU PHOTONICS UK LIMITED: 2 Howard Court, 10 Tewin Road, Welwyn Garden City, Hertfordshire, AL7 1BW, UK, Telephone: (44)1707-294888, Fax: (44)1707-325777 E-mail: [info@hamamatsu.co.uk](mailto:info@hamamatsu.co.uk)

North Europe: HAMAMATSU PHOTONICS NORDEN AB: Torshamnsgatan 35, 16440 Kista, Sweden, Telephone: (46)8-509-031-00, Fax: (46)8-509-031-01 E-mail: [info@hamamatsu.se](mailto:info@hamamatsu.se)

Italy: HAMAMATSU PHOTONICS ITALIA S.R.L.: Strada della Moia, 1 int. 6 20044 Arese (Milano), Italy, Telephone: (39)02-93 58 17 33, Fax: (39)02-93 58 17 41 E-mail: [info@hamamatsu.it](mailto:info@hamamatsu.it)

China: HAMAMATSU PHOTONICS (CHINA) CO., LTD.: 1201, Tower B, Jiaming Center, 27 Dongsanhuan Bellu, Chaoyang District, 100020 Beijing, P.R. China, Telephone: (86)10-6586-6006, Fax: (86)10-6586-2866 E-mail: [hpc@hamamatsu.com.cn](mailto:hpc@hamamatsu.com.cn)

Taiwan: HAMAMATSU PHOTONICS TAIWAN CO., LTD.: 13F-1, No.101, Section 2, Gongdao 5th Road, East Dist., Hsinchu City, 300046, Taiwan(R.O.C) Telephone: (886)3-659-0080, Fax: (886)3-659-0081 E-mail: [info@hamamatsu.com.tw](mailto:info@hamamatsu.com.tw)